



Making advancements in rotator cuff
repair and treatment

 **smith&nephew**
MAGNUM[®] X
Stainless Steel Knotless Implant




MAGNUM[◇] X

Stainless steel knotless implant

- Independent suture tensioning gives the surgeon directional control of the repair tension
- Bone lock is achieved before suture lock, providing the ability to confirm final implant location before final tension is locked
- Internal suture lock eliminates reliance on bone quality for suture lock, maintaining final tension over time
- Knotless design locks securely under the cortex, reducing the possibility of gap formation

MAGNUM X implant specifications

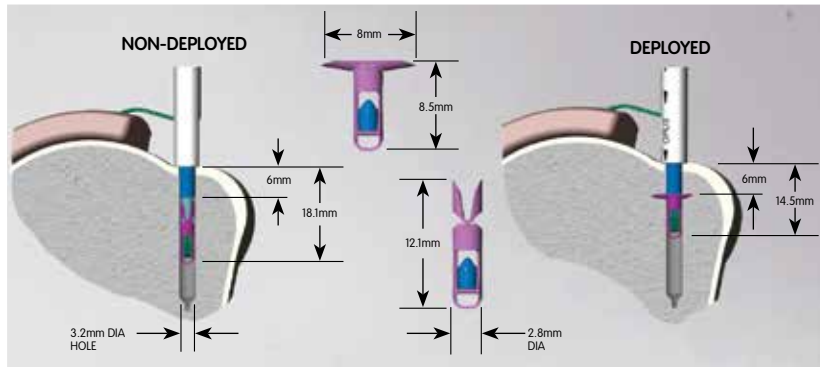

SMARTSTITCH[◇] suturing system

The SMARTSTITCH system is an integrated grasper and suturing mechanism that allows the surgeon to grasp the cuff tissue and quickly place an Incline Mattress, or side-to-side stitch. The SMARTSTITCH system consists of a reusable SMARTSTITCH handle and a single-use PERFECTPASSER[◇] insert.



Ordering information

MAGNUM X stainless steel knotless implant

Reference #	Description
OM-1503	MAGNUM X 2.8MM KNOTLESS IMPLANT
OM-3001	MAGNUM ² INSTRUMENT SET: INCLUDES OM-8000, OM-9010, OM-9015, OM-9020, OM-9330
OM-9136	MAGNUM X PLUS 2.8MM KNOTLESS IMPLANT SET WITH PERFECTPASSER MAGNUMWIRE [◇] SUTURE CARTRIDGE, WHITE
OM-9137	MAGNUM X PLUS 2.8MM KNOTLESS IMPLANT SET WITH PERFECTPASSER MAGNUMWIRE SUTURE CARTRIDGE, BLUE CO-BRAID
OM-9010	3.2MM SPADE DRILL
OM-9015	DRILL GUIDE AND PATHFINDER [◇] SET: COMPATIBLE WITH OM-9010, OM-9020, OM-9410 AND OM-9420
OM-9020	3.3MM BONE PUNCH

ArthroCare Corporation
7000 West William Cannon Drive
Austin, TX 78735
USA

www.smith-nephew.com

Order Entry: 1-800-343-5717
Order Entry Fax: 1-888-994-2782

© 2015 Smith & Nephew, Inc.
*Trademark of Smith & Nephew. Reg. US Pat. & TM Office.

P/N A1075 Rev. F 01/15