


 smith&nephew

**XL
Microfracture
Pick**

90°



90 Degrees of Precision

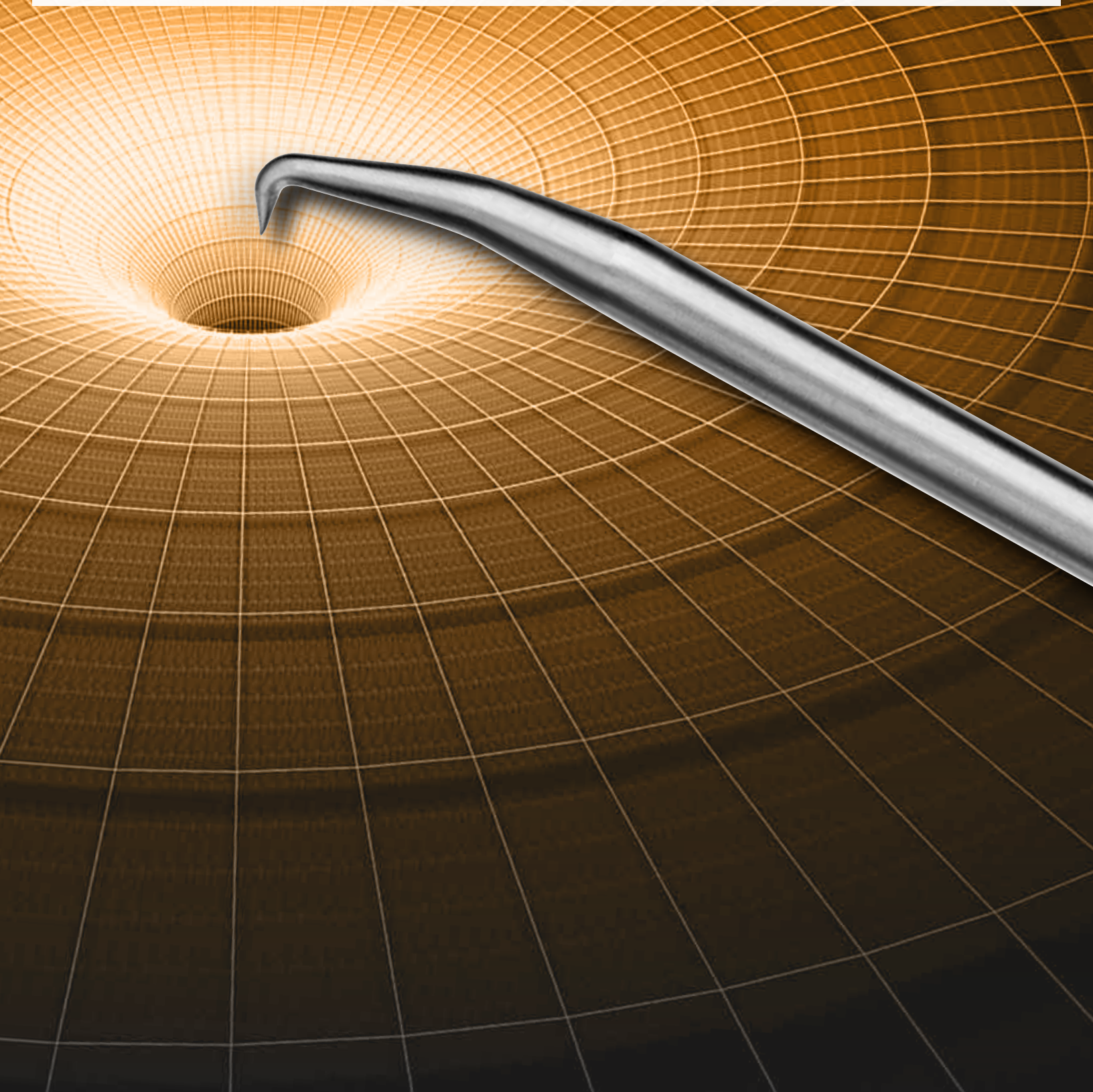

Applies force where you need it.

Designed for treatment of chondral lesions in the hip joint, the XL MICROFRACTURE PICK has a strategically-placed strike point that enables the mallet force to be transferred in the direction of the 90-degree pick. This allows the creation of a precisely located hole, while minimizing skiving and iatrogenic damage.

Optimally-positioned strike point delivers force in the direction of the mallet strike, minimizing skiving and iatrogenic damage.

High strength Delrin® resin handle provides a secure and ergonomic grip.

Tapered shaft with 15.2" working length enables access through a standard arthroscopy portal, using a slotted cannula.

Sharp needle-point is angled at 90 degrees to allow force to be delivered at the site of the chondral lesion.



Image courtesy of Marc J. Philippon, MD



Ordering Information

XL Microfracture Pick

Reference #	Description
72203338	XL Microfracture Pick, 90 Degree

Smith & Nephew, Inc.
150 Minuteman Road
Andover, MA 01810
USA

www.smith-nephew.com
T +1 978 749 1000
US Customer Service: +1 800 343 5717
International Customer Service: +1 978 749 1140

®Trademark of Smith & Nephew. Registered US Patent and Trademark Office. ©2013 Smith & Nephew. All rights reserved. Printed in USA. 04/13 10601093 Rev. A